


Operator's Manual

Chart Projector



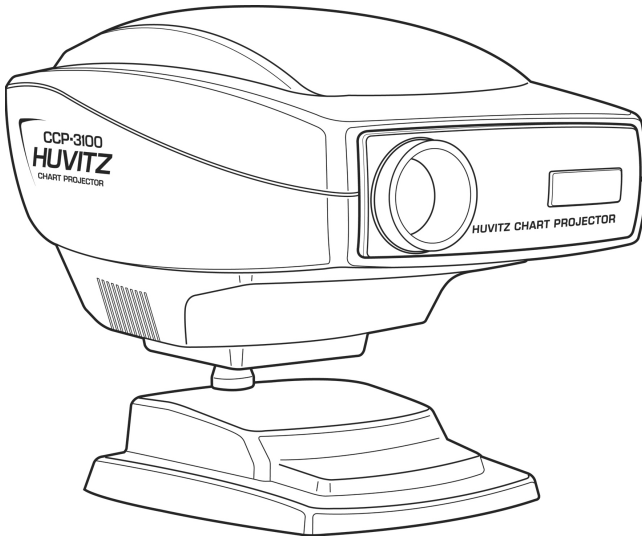

IMPORTANT NOTICE

This product may malfunction due to electromagnetic waves caused by portable personal telephones, transceivers, radio-controlled toys, etc. Be sure to avoid having objects such as, which affect this product, brought near the product.

The information in this publication has been carefully checked and is believed to be entirely accurate at the time of publication. HUVITZ assumes no responsibility, however, for possible errors or omissions, or for any consequences resulting from the use of the information contained herein.

HUVITZ reserves the right to make changes in its products or product specifications at any time and without prior notice, and is not required to update this documentation to reflect such changes.

All rights are reserved.

Under copyright laws, this manual may not be copied, in whole or in part, without the prior written consent of HUVITZ Co., Ltd.

Document Number: 90000MCP001-A, published in 2006/07/20

CONTENTS

1. Introduction	5
1.1. Outline of product	5
1.2. Classifications	5
1.3. Certification	5
2. Safety Information	6
2.1. Introduction.....	6
2.2. Safety Symbols.....	7
2.3. Environmental factors	8
2.4. Safety Precautions	9
3. Features.....	11
4. Notes for Using the Instrument	12
5. Configurations.....	13
5.1. Main Body	13
5.2. Main Body cont.....	14
5.3. Accessories	15
6. Installation..	16
6.1. Position of Chart Projector.....	16
6.2. How to obtain correct Letter Size and focus (Zooming Type only).....	17

4	CCP-3100 Operator's Manual	
7.	Input and Use of Program	18
7.1.	How to input user's program	18
7.2.	How to use your stored programs.	19
8.	Types of Chart (Type C is Standard)	20
9.	Maintenance and Repair	21
9.1.	Lamp Replacement.....	21
9.2.	Fuse Replacement.....	22
9.3	The battery of Remote Controller Replacement	23
9.4	Cleaning.....	24
9.5	Disposal	24
10.	Specifications	25
11.	Components List	26
11.1.	Standard accessories	26
11.2.	Optional accessories	26
12.	Service Information	27

1. Introduction

1.1. Outline of product

HUVITZ CCP-3100 is an automatic chart projector which provides 41 different charts including red/green and polarized filters for the eye optometry. Users can select the desired charts instantaneously by cordless remote controller. Charts move quickly and quietly. The customized two programs will help you to perform more convenient and effective examination. You can also execute the binocular balance tests, stereo tests, fixation disparity tests, fusion tests and aniseikonia tests with CCP-3100. This automatic chart project is compatible with HUVITZ digital refractor, CDR-3100.

Equipment intended for connection to signal input and signal output, shall with IEC60601-1, or UL60601-1 harmonized national standard or the combination. If, in doubt, contact qualified technician for help.

1.2. Classifications

- Class I equipment (CCP-3100)
- No applied part
- Ordinary equipment (IPX 0)
- Equipment not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide
- Continuous operation equipment

1.3. Certification



UL60601-1
CAN/CSA C22.2
NO.601.1



**MEDICAL EQUIPMENT
WITH RESPECT TO ELECTRIC SHOCK
FIRE AND MECHANICAL HAZARDS ONLY
IN ACCORDANCE WITH UL 60601-1, AND
CAN/CSA C22.2 NO.601.1**

2.Safety Information

2.1 Introduction

Safety is everyone's responsibility. The safe use of this equipment is largely dependent upon the installer, user, operator, and maintainer. It is imperative that personnel study and become familiar with this entire manual before attempting to install, use, clean, service or adjust this equipment and any associated accessories. It is paramount that the instructions contained in this manual are fully understood and followed to enhance safety to the patient and the user/operator. It is for this reason that the following safety notices have been placed appropriately within the text of this manual to highlight safety related information or information requiring special emphasis. All users, operators, and maintainers must be familiar with and pay particular attention to all Warnings and Cautions incorporated herein.

WARNING

“Warning ” indicates the presence of a hazard that could result in severe personal injury, death or substantial property damage if ignored.

NOTE

“Note” describes information for the installation, operation, or maintenance of which is important but hazard related if ignored.






CAUTION

“Caution” indicates the presence of a hazard that could result in minor injury, or property damaged if ignored.

2.2 Safety Symbols




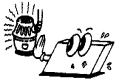




The International Electrotechnical Commission (IEC) has established a set of symbols for medical electronic equipment, which classify a connection or warn of any potential hazards. The classifications and symbols are shown below.

Save these instructions

	<p>I and O on the power switch represent ON and OFF respectively.</p>
	<p>This symbol identifies a safety note. Ensure you understand the function of this control before using it. Control function is described in the appropriate User's or Service Manual.</p>
	<p>Hot surface.</p>
	<p>Identifies the point where the system safety ground is fastened to the chassis. Protective earth connected to conductive parts of Class I equipment for safety purposes.</p>
	<p><u>Disposal of your old appliance</u></p> <ol style="list-style-type: none"> 1. When this crossed-out wheeled bin symbol is attached to a product it means the product is covered by the European Directive 2002/96/EC. 2. All electrical and electronic products should be disposed of separately from the municipal waste stream via designated collection facilities appointed by the government or the local authorities. 3. The correct disposal of your old appliance will help prevent potential negative consequences for the environment and human health. 4. For more detailed information about disposal of your old appliance, please contact your city office, waste disposal service or the shop where you purchased the product.

2.3 Environmental factors

Avoid the following environments for operation or storage:

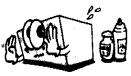	<p>Where the equipment is exposed to water vapor. Don't operate an equipment with a wet hand.</p>
	<p>Where the equipment is exposed to direct sunlight.</p>
	<p>Where the temperature changes extremely. Normal operating temperature range is from 10°C to 40°C, Humidity is from 30% to 75%.</p>
	<p>Where it is near the heat equipment.</p>
	<p>Where the humidity is extremely high or there is a ventilation problem.</p>
	<p>Where the equipment is subject to excessive shocks or vibrations.</p>
	<p>Where equipment is exposed to chemical material or explosive gas.</p>
	<p>Be careful not to be inserted dust, especially, metal</p>

Avoid places where the ambient temperature falls below 10°C or exceeds 40°C for normal operation, or below -10°C or exceeds 40°C (14°F-104°F) for transportation and storage. Humidity should be maintained between 30 and 75% for normal operation, transportation and storage and atmospheric pressure (700~1060hpa). Avoid environments where the equipment is exposed to excessive shocks or vibrations.

2.4 Safety Precautions

This equipment has been developed and tested according to national as well as international safety standards. This guarantees a very high degree of safety for this device. The correct handling of this equipment is imperative for its safe operation. Therefore, please read carefully all instructions before switching on this device. For more detailed information, please contact our Customer Service Department or one of our authorized representatives.

1. This equipment must not be used (a) in an area that is in danger of explosions and (b) in the presence of flammable, explosive, or volatile solvents such as alcohol, benzene or similar chemicals.
2. Do not put or use this device in humid rooms. Humidity should be maintained between 30 and 75% for normal operation. Do not expose the device to water splashes, dripping water, or sprayed water. Do not place containers containing fluids, liquids, or gases on top of any electrical equipment or devices.
3. The equipment must be operated only by, or under direct supervision of a properly trained and qualified person.
4. Customer maintenance of this equipment may only be performed as stated in the User's Manual. Any additional maintenance may only be performed by HUVITZ service technicians or other authorized persons.

10 CCP-3100 Operator's Manual -----

5. The manufacturer is only responsible for effects on safety, reliability, and performance of this equipment when the following requirements are fulfilled: (1) The electrical installation in the respective room corresponds to the specifications stated in this manual and (2) This equipment is used, operated, and maintained according to this manual.
6. The manufacturer is not liable for damage caused by unauthorized tampering with the device(s). Such tampering will forfeit any rights to claim under warranty.
7. Only persons who have undergone proper training and instructions are authorized to install, use, operate, and maintain this equipment.
8. Keep the User's Manual in a place easily accessible at all times for persons operating and maintaining the equipment.
9. Please do not pull on any cable. Always hold on to the plug when disconnecting cables.
10. Before every operation, visually check the equipment for exterior mechanical damage(s) and for proper function.
11. Do not cover any ventilation grids or slits.
12. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to order devices in the vicinity. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures :
 - Reorient or relocate the receiving device.
 - Increase the separation between the equipment.
 - Connect the equipment into an outlet on a circuit different from that to which the other device are connected.
 - Consult the manufacturer or field service technician for help.

3.Features

Bright (halogen), high quality images

Quiet, instantaneous chart selection

41 different charts including: red/green and polarizing filters, binocular balance tests, stereo tests, fixation disparity tests, fusion tests and Aniseikonia tests.

34 Masks (vertical line, horizontal line, and single character)

2.5 ~ 8 m (8 ~ 26 feet) projection range

Cordless Remote Control

User Programmable (2 programs up to 30 charts each)

Easily networked with HUVITZ's Digital Refractor

4. Notes for Using the Instrument

Do not hit or drop the instrument. The instrument may be damaged if it receives a strong impact. The impact can damage the function of this instrument. Handle it with care.

Exposure to direct sunlight or very bright indoor lights can influence the results of measurements.

If you want to connect this with other equipment, consult the dealer.

Sudden heating of the room in cold areas will cause condensation of vapor on the protective glass in the measurement window and on optical parts inside the instrument. In this case, wait until the condensation disappears before performing measurements.

Keep the objective glass of the examinee side clean. If it is stained with other substances or a fingerprint, it may cause an error or inaccurate measurements.

Disconnect the power supply and consult the dealer when there is smoke, strange odors, or noise while working.

Don't use organic solution such as alcohol, thinner, benzene, etc. to clean the surface of this instrument. It may damage the instrument.


When moving the CCP-3100, always make sure the power supply is off, and then lift the bottom of the unit with both hands.

If you leave the CCP-3100 without using it for long periods, disconnect the power supply and protect the unit with dust cover.

5. Configurations

5.1 Main Body

The CCP-3100 Chart Projector is mainly composed of the body, stand, screen and remote control.

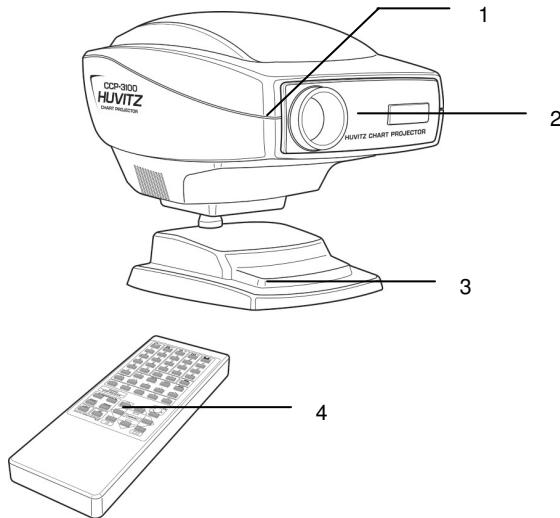


[Fig. 1 Front View with Remote Control]


- 1. Projection lens
- 3. Stand (table mount)

- 2. IR receiving window
- 4. Wireless remote control

5.2 Main Body cont.

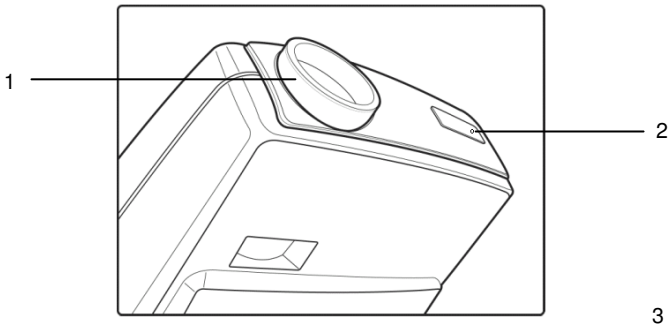


[Fig. 2 Bottom View]

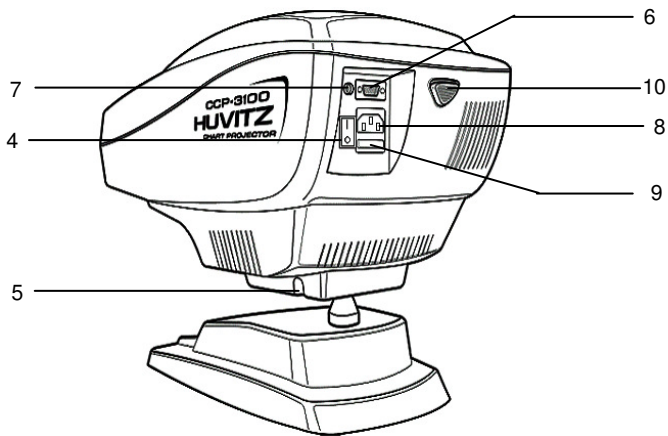

[Fig. 3 Rear View]

-
- | | | |
|----------------------------------|-----------------------|------------------------|
| 1. Focusing Wheel | 2. Power ON/OFF LED | 3. 110/220 Volt Switch |
| 4. Power On/Off Switch | 5. Tilting Adjustment | 6. RS-232 Comm.Port |
| 7. CAN (Controller Area Network) | 8. AC Power Inlet | 9. Fuse Holder |
| 10. Cover Release | | |

5.3 Accessories

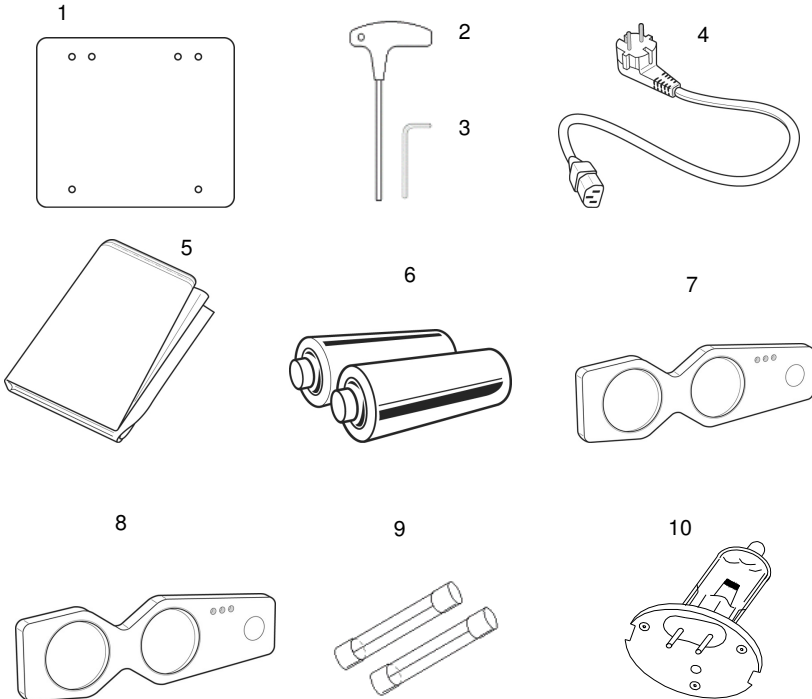


[Fig.4 Accessories]

-
- | | | |
|--|-----------------------|------------------------------------|
| 1. Screen | 2. Allen wrench (5mm) | 3. Allen wrench (2.5mm) |
| 4. Power cable | 5. Dust Cover | 6. Batteries (2 size AA(LR6) 1.5V) |
| 7. Polaroid glasses (Optional) | | |
| 8. Red/Green glasses (Optional) | | |
| 9. Fuse | | 10. Lamp (6V20W) |
| 9.1 Primary Fuse (2 size 5X20mm 250V T0.63AL) | | (6V30W)--Optional |
| 9.2 Secondary Fuse (F1 : 1 size 5X20mm 250V T1AL) | | |

6.Installation

6.1 Position of Chart Projector

The single tube projector must be installed at the same distance from the screen as the patient.

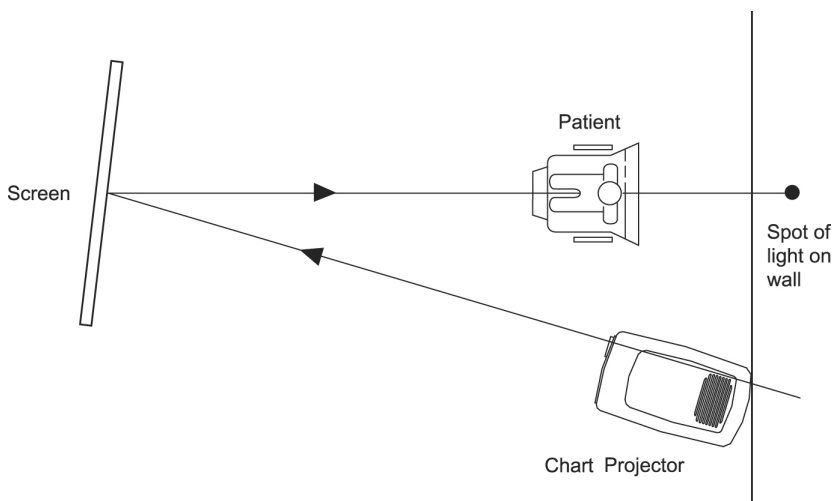

[Fig. 6 Installation without mirrors]

1. Install chart projector at desired position within 2.5~8m (8~26 feet).
2. Set the projection angle of chart processor based on Screen.
3. Project 0.05(20/400, 6/120) "C" chart onto the screen.
4. Adjust it for a bright clear image using the focusing wheel.

6.2 How to obtain correct Letter Size and focus (Zooming Type only)

1. Install chart projector at desired position of 2.5~8m (8~26 feet)
2. Measure the distance from patient to screen
3. Refer to the test scale shown in the Appendix A
4. Have someone hold the test scale at the screen
5. Project the 20/400(6/120, 0.05) alphabet "C"
6. Adjust the position of the projector forward or backward and turn the focusing wheel until "C" is brought in sharp focus within the respective distance scale

7. Input and Use of Program

7.1 How to input user's program

- Step 1. Press the "LIGHT" button **once**.
Results: The light goes out.
- Step 2. Press the "PGM 1" button **three** times.
Results: The light will come on and the chart projector will display the first chart of the default program.
- Step 3. Select the first chart for your new program.
Press the "Program \blacksquare " button **once**.
- Step 4. Select the second chart for your new program.
Press the "Program \blacksquare " button **once**.
- Step 5. Repeat steps 3 & 4 as desired for a total of thirty charts per program.
Note: If you want to check or modify the previously programmed chart, press the "Program \blacksquare " button.
- Step 6. Press the "PGM 1" button **one** time to save your new program.
Caution : If you do not press the PGM1 button, your program will not be saved.

Program No. 2


Repeat Steps 1 through 6 substituting "PGM 2" for "PGM 1".

7.2 How to use your stored programs.


- Step 1. Press " PGM 1 " or "PGM 2" to select the program of your choice.
- Step 2. Use the "Program ▢" or "Program ▢" to scroll through your stored program.

8.Types of Chart (Type C is Standard)


A TYPE

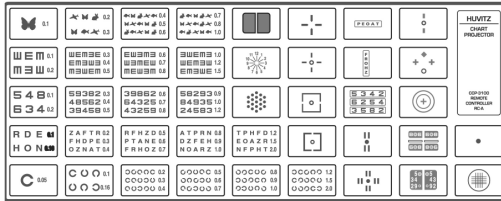

C TYPE

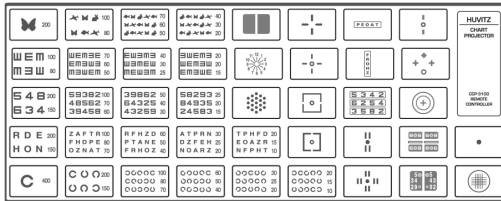

D TYPE

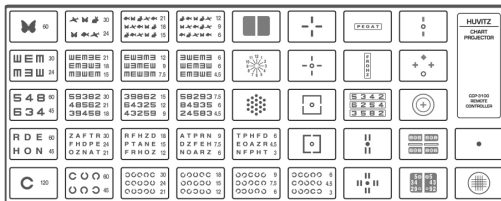

[Fig. 7 Types of Chart]

9.Maintenance and Repair

Maintenance : There is no periodic or routine user maintenance required.

9.1 Lamp Replacement


1. Turn off power and disconnect the power cord.
2. Allow the projector to cool down for about 5 minutes before replacing the lamp.
3. Using the 2.5mm Allen wrench remove the two screws from the bottom front of the projector cover.
4. Press the button at the rear of the projector and remove the top cover.
5. Remove the lamp attaching screws, hold socket and remove the lamp.
6. Hold the socket and the lamp bracket and remove the lamp from the socket.
7. Replace the lamp. (Fig.8 Lamp assembly)
8. Connect the power cord and check the lamp.
9. Reinstall the top cover.



WARNING


Do not touch the lamp while it is hot. (Risk of burns)

Do not touch the projector with wet hands. (Risk of electric shock)

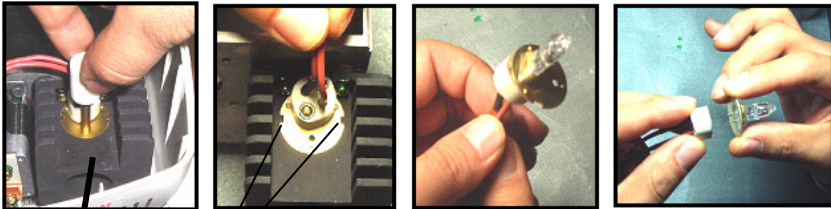

[Fig. 8 Lamp assembly]


**Install the lamp securely as shown in the Fig.8
Make sure that you install it in the correct direction.**

 **WARNING HOT**
Turn off power and disconnect the power cord. Allow the projector to cool down for about 5 minutes before replacing the lamp because the heated lamp may burn your fingers.

9.2 Fuse Replacement

1. Turn off power, and disconnect the power cord.
2. Pull out the plastic cover at the lower part of power socket, to replace the fuse.

 **WARNING**
Do not touch the projector with wet hands. (Risk of electric shock)

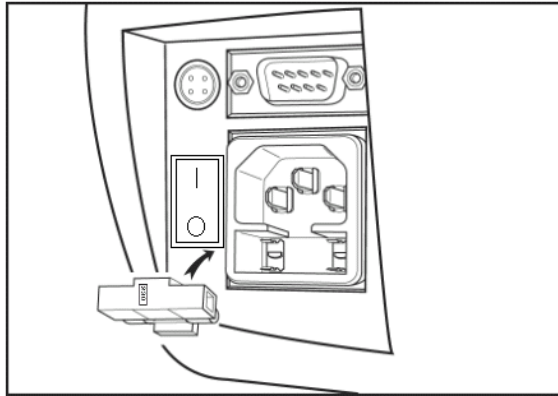

[Fig. 9 Fuse assembly]

9.3 Replacement the Battery of Remote Controller

1. Remove the battery cover of the remote controller.
2. Replace the batteries and the battery cover to the remote controller.

NOTE : Be sure to observe the polarity marks (plus + & minus -) inside the battery compartment.

Remove these batteries if equipment is not likely to be used for some time. Follow local governing ordinances and recycling plans regarding disposal or recycling of device components.



Battery(2)
LR6(AA)1.5V

9.4 Cleaning

There are no cleaning requirements other than regular office housekeeping, such as dusting.

1. Wipe the cover and the screen with a dry soft cloth when they get dirty.
2. If for some reason the projector should become soiled, wipe it clean with the damp lint-free cloth and mild detergent.
3. When the lens is dirty, wipe it with a dry soft cloth.

9.5 Disposal

NOTE

To dispose the instrument, accessories, and components, follow local governing ordinances and recycling plans regarding disposal or recycling of instrument or device components. Especially a lithium battery may pollute the environment if the instrument or a lithium battery is abandoned.

When disposing packing materials, sort them by the materials and follow local governing ordinances and recycling plans.

10.Specifications

Charts	41charts -34 mask, Red/green& Polarization Filters
Distance of projection	2.5~8m(Optional: 8m screen)
Chart Rotation Speed	Average 0.15 sec
Chart magnification	30X (When 5m distance)
Power saving	Automatic Lamp Off (10 min)
Program	2 programs with a maximum of 30 charts each
Tilt angle	+/- 15 degrees
Power Supply	100-120V~, 200-240~, 50/60 Hz, 0.6A
Lamp	6V/20W , 6V/30W(Optional)
Fuse	250V, T0.63AL
Dimension	367X219X273mm
Weight	6.6kg (Body: 5.8kg + Stand: 0.8kg)

11.Components List

11.1 Standard accessories

1. Remote control -1
2. Screen (400X500mm) -1
3. Table stand -1
4. Dust cover -1
5. Power cable –1

For use of equipment in rated voltage less than 125Vac,minimum 6A,Type SJT or SVT , 18/3AWG,10A, max 3.0m long : One end with Hospital Grade Type, NEMA 5-15P Other end with appliance coupler.

For use of equipment in rated voltage less than 250Vac,minimum 6A,Type SJT or SVT , 18/3AWG,10A, max 3.0m long : One end terminated with blade attachment plug(HAR) Type, NEMA 6-15P.
6. Allen wrenches (2.5mm and 5mm) -2
7. Fuse – 3
8. Lamp – 1
9. Batteries -2

11.2 Optional accessories

1. Red/green glasses
2. Polarization glasses

12. Service Information

How to contact service: If there are any problems with the equipment, please follow the steps below:

- First of all, refer to the 9. Maintenance and Repair sections according to the problem that you are encountered. And then follow the suggested sequences.
- If the problem persists, please contact the local distributor in your province or country at first.
- Before calling to the local distributor, you'd better check these information such as Model and Serial Numbers. To do so, fill up the following table as soon as you purchase our product. You can look up these information at any time. The serial number is found on the back of this unit. The serial number is unique to this unit. You should retain this manual as a permanent record of your purchase. Please retain your purchase receipt as your proof of purchase.

Date of Purchase:

Dealer's Name:

Dealer Address:

Dealer Phone No.:

Model No.:

Serial No.:


■ If you can't contact with your local distributor, you can directly get in touch with the service department of the HUVITZ using the phone number and the address written in the below table.

How to Contact HUVITZ Co., Ltd


Write us at:	Tel: +82-31-442-8868
HUVITZ Co., Ltd.	Fax: +82-31-442-8619
Huvitz B/D, 689-3 Geumjeong-dong	URL: http://www.huvitz.com
Gunpo-si Gyeonggi-do, South Korea	e-mail: ccp-support@huvitz.com
435-862	


■ **Appendix A**



Test Scale Hold this test scale at the Screen, and according to the refracting distance from the patient's eye to the screen such as 3m, 4m, 5m or 6m position the projector backward and forward to fit the projector Alphabet 400 chart "C" with this Scale frame in sharp focus

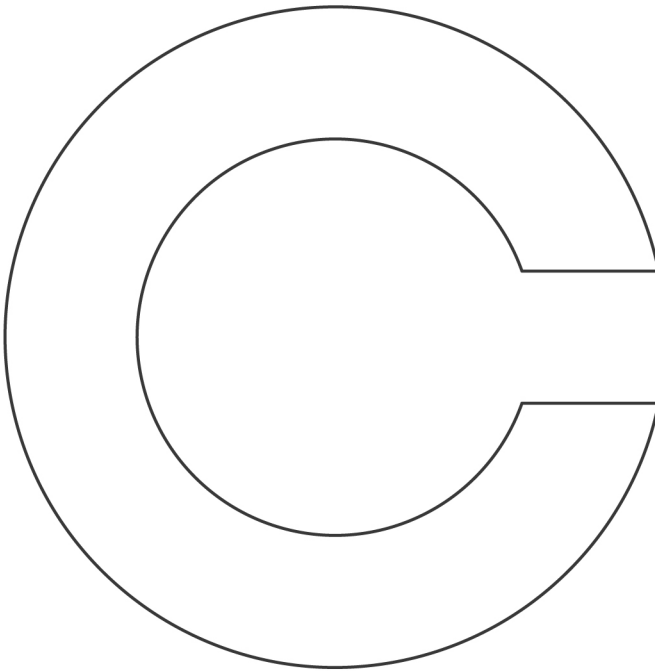

3 m

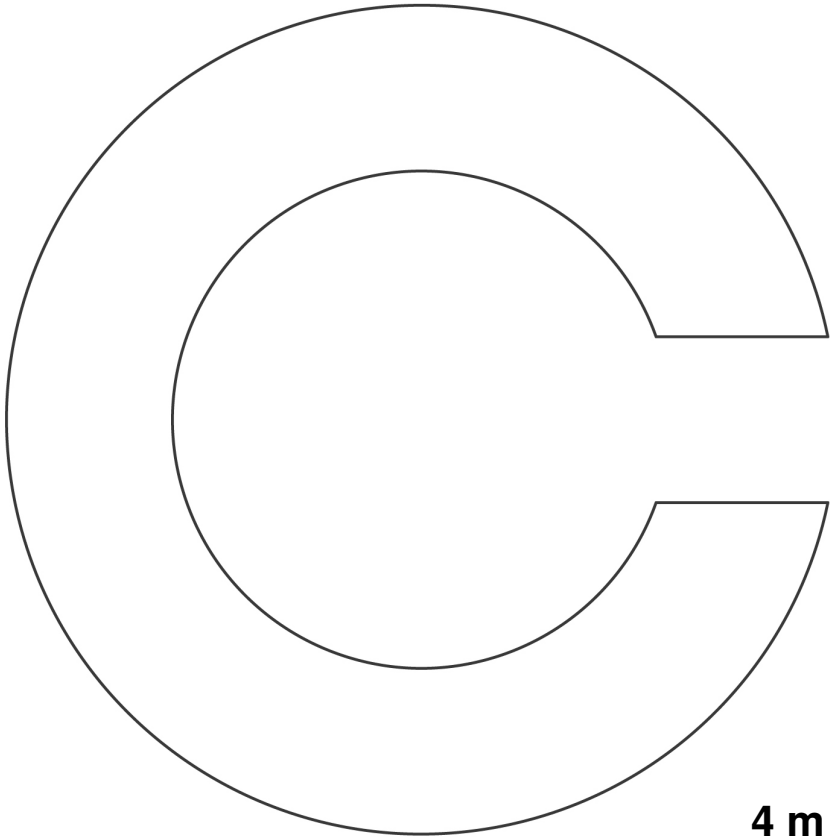

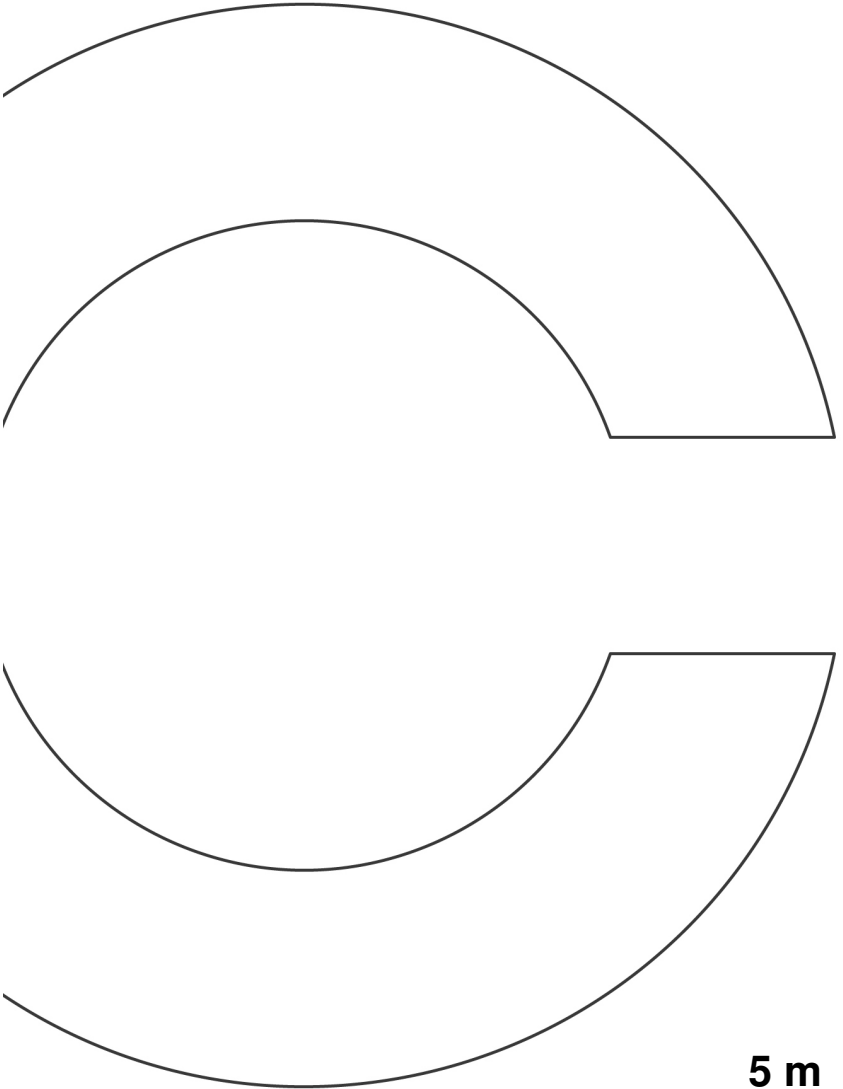

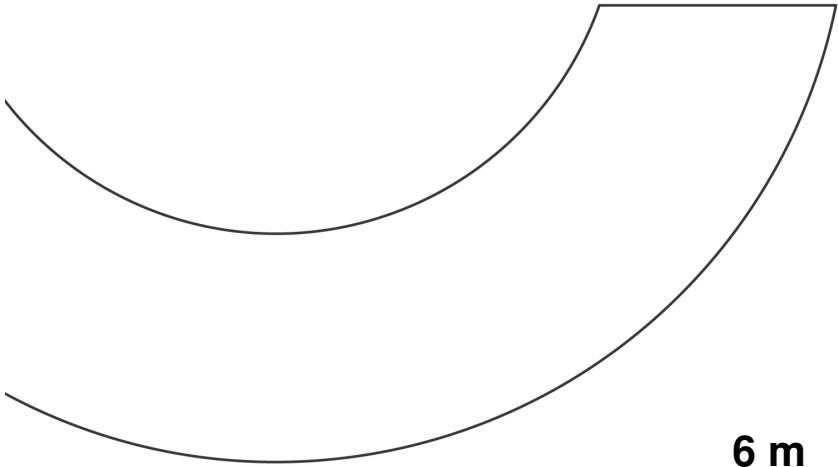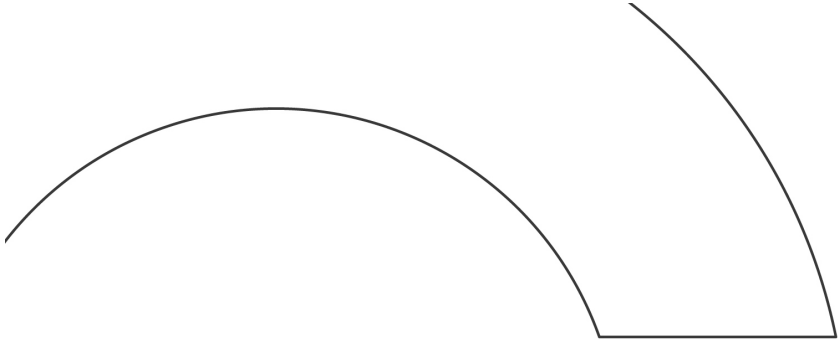

6 m